EKULF

CORDLESS PREMIUM POWERFLOSSER

INSTRUCTION MANUAL

**EXTRA BIG
WATER TANK**

**WATER PROOF
DESIGN**

**OPERATION
MODES**

**0,65 mm
WATER FLOW**

**SMART POWER
DISPLAY**

CORDLESS PREMIUM POWERFLOWER

FLOSS WITH THE POWER OF WATER

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IMPORTANT SAFEGUARDS

All information in this section is provided to instruct you in the correct and safe operation of this device so that you can avoid injuries to yourself and others as well as damage to property.

Before operating this product, please read these instructions completely and save this manual for future use.

READ ALL INSTRUCTIONS BEFORE USING!

DANGERS


- Toddlers, infants, and people who are unable to operate should not be allowed to use this appliance. Please keep it out of reach of children.
- Never alter the appliance. Also, do not try to disassemble or repair it.
- Do not place or store the product where it can fall or be pulled into a tub or sink.
- Do not reach for the product that has fallen into water. Unplug immediately.
- Do not use the product while bathing.

WARNINGS

This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory, or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

- Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- The appliance is only to be used with the power supply unit provided with the appliance.
- The battery must be removed from the appliance before it is scrapped.
- The appliance must be disconnected from the supply mains when removing the battery.
- The battery is to be disposed of safely.
- Regarding how to remove the battery, please see "Disposal" on page 27.
- Do not fill the tank with water warmer than 40 °C.
- If the external flexible cable or cord of this transformer is damaged, it shall be replaced by the manufacturer, his service agent, or a similarly qualified person in order to avoid a hazard.
- Keep the cord away from heated surfaces.
- Please do not use it if you lose the nozzle tip.
- Do not submerge the device in water.
- Do not use for any purpose other than oral cleaning.
- Keep the nozzle tip away from small children. The nozzle tip is small enough to be swallowed and can become lodged in the throat.
- Certain types of mouthwash may damage the equipment by causing the housing or water reservoir to crack. Therefore, do not use mouthwash in the water reservoir.
- People who experience tooth or gum pain may not be able to use this device and should consult their dentist.
- Do not press the nozzle tip too firmly against your teeth or gums. Doing so may injure your gums.

SAVE THESE INSTRUCTIONS

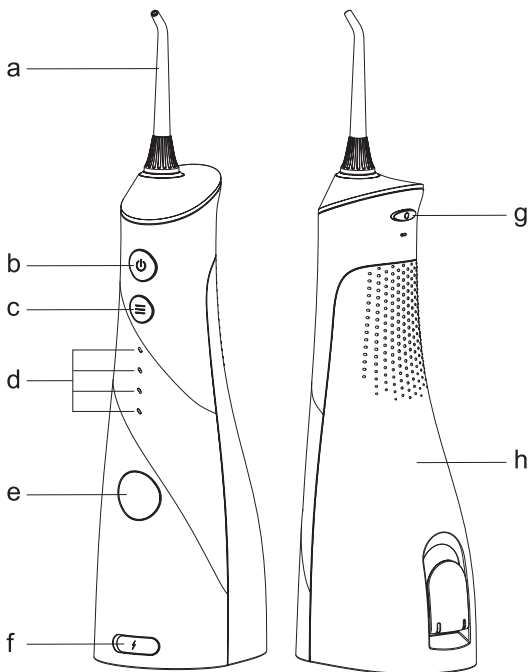CONFIGURATION

- a. Nozzle tip
- b. On/Off button
- c. Mode button
- d. Mode indicator

- e. Screen display
- f. Charging port
- g. Nozzle tip release button
- h. Water tank

Note: There may be some liquid remaining inside the main appliance before use (tank, water suction tube, etc). This is the remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.

HOW TO USE

Charging the device (Fig. 1)

1. Take the charger base cover down and insert the appliance plug securely into the appliance socket (make sure no liquid remains inside the socket before inserting the plug).
2. Insert the charging adapter securely into the power outlet.
3. When the device is charging, the battery indicator light on the display screen lights up and flashes to show real-time battery level. Once fully charged, the light remains steady for 10 seconds before turning off.
4. Replace the charger base cover after charging is finished.

Note: If the power indicator light on the display screen starts blinking, it indicates that the battery needs to be charged. A full charge will take about 4 hours.

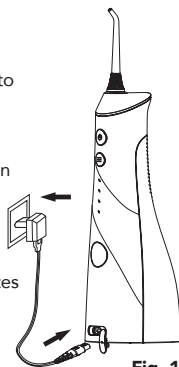
Fig. 1

Use the Oral Irrigator

1. Attach the nozzle tip onto the unit. (Fig. 2)
2. Open the water tank cap to fill the tank with water at a temperature not exceeding 40 °C, then close the cap. (Fig. 3)
3. Or, slide and remove the water tank, fill the tank with water not exceeding 40 °C, then attach the tank to the unit. (Fig. 4)

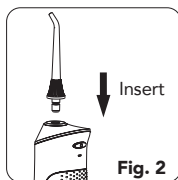
Fig. 2

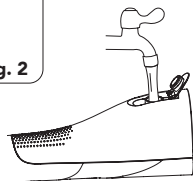
Fig. 3

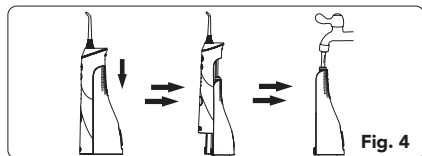
Fig. 4

4. Press the "Mode" button to select the desired mode. During irrigation, you can change the working mode:

Clean: for normal cleaning

Soft: for gentle cleaning of sensitive gums

Manual: Press the button to start cleaning, release the button to stop

Customize: Press and hold the "Mode" button for 1.5 seconds to enter custom mode. Each press increases the level by one (cycling from 1 to 8 and back to 1). The current level will be locked if no further presses are made within 5 seconds.

5. Hold the unit upright with the nozzle tip pointed toward your teeth.
6. Keep your mouth slightly open so that the water can flow out.
7. Press the "On/Off" button to start irrigating:
 - a. Direct the water jet flow at a 90-degree angle to the teeth and gums.
 - b. Guide the water jet slowly along the teeth.
 - c. Direct the water jet along the gum line with the nozzle tip at a 90-degree angle to the gum tissue.
 - d. Clean between the teeth and around any braces, crowns, or bridgework thoroughly. (Fig. 5)

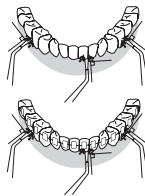
Fig. 5

8. Please turn off the irrigator before replacing the nozzle tip, then press the nozzle tip release button to replace it with another nozzle tip. (Fig. 6)
9. Once irrigation is complete, press the "On/Off" button again to turn off the device.

After Using the Appliance (Fig. 7)

1. Remove the water tank and drain the water.
2. Press the "On/Off" button to turn on the unit and drain remaining water inside.
3. Press the "On/Off" button again to turn off the unit and close the water tank cap.
4. Wipe the device with a dry cloth.

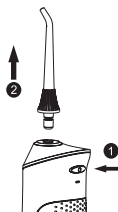
Fig. 6

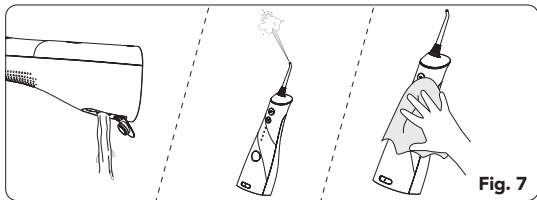
Fig. 7

Note: To prevent bacterial growth, ensure no water remains inside the unit. If you're not going to use this appliance for a long time, be sure to wipe and dry the appliance before storing it.

MAINTENANCE

Use only water or a neutral cleaning agent to clean the device. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), as they may damage the device. Do not wash using hot water over 40 °C. Keep the device away from high temperatures and avoid direct sunlight.

Main unit

1. Remove dirt by wiping with a cloth and liquid soap, then dry with a clean cloth.
2. Do not submerge the main unit in water.
3. Dry thoroughly after cleaning.

Nozzle tip

1. Rinse with water and wipe with a soft cloth.
2. Replace the nozzle tip every 6 months.

Water suction tube

1. Rinse with water and wipe with a soft cloth.
2. Do not bend, pull, or twist the tube.

Water tank

1. Rinse with water and wipe with a soft cloth.
2. Dry the tank if not using it for an extended period.

DISPOSAL

The rechargeable battery inside your appliance is easy to remove for recycling. At the end of the product's life and prior to disposal, remove the battery from the handle by following the instructions below. Be sure to protect your eyes, hands, fingers, and the work surface.

1. Slide the water tank off the device and unscrew the two screws beside the suction tube. (Fig. 8)
2. Remove the screw from the bottom of the handle. (Fig. 9)
3. Remove the fixing bracket and PCBA, then melt the tin and take the battery out. (Fig. 10)
4. Recycle or dispose of the battery properly. Ask your local waste management office if they accept the remaining components for recycling. (Fig. 11)

If you have recycling questions, please contact your local waste management office.

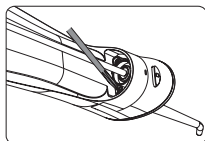
Fig. 8

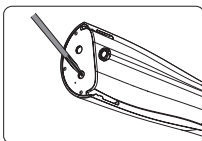
Fig. 9

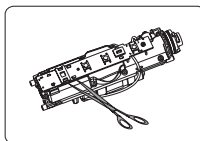
Fig. 10


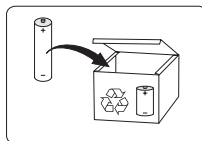
Fig. 11

TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
It does not work.	The product is newly purchased or unused for more than 3 months.	Charge the unit for at least 4 hours.
It only works a few minutes, even after charging.	The battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the unit for at least 4 hours.
Water pressure is very low.	Dirt is stuck in the nozzle tip.	Replace the nozzle tip.
	The nozzle tip is deformed.	
	Dirt is stuck in the filter.	Clean the filter.
Water does not come out.	The tank is empty.	Fill the tank.
	The appliance is tilted too much.	Keep appliance upright during use.
	The product is newly purchased or unused for a long time.	Remove the nozzle tip from the appliance and flush the nozzle hole.

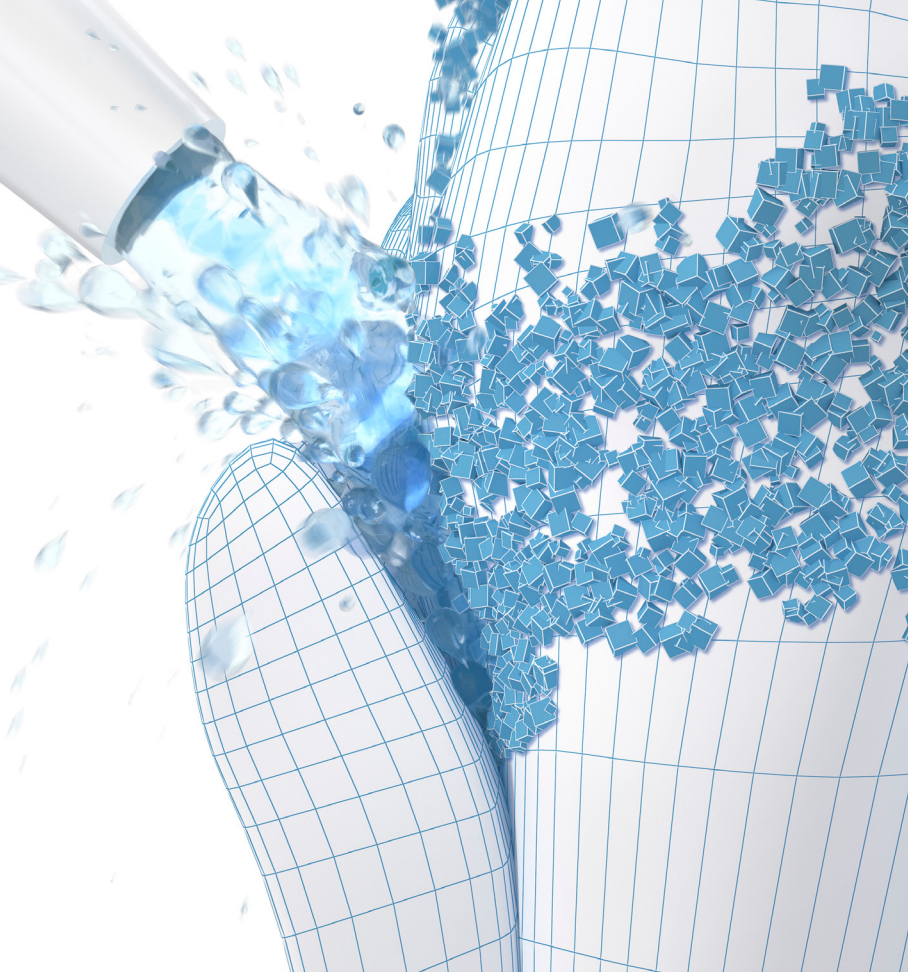CORDLESS PREMIUM
POWERFLOSSER

FLOSS WITH THE POWER OF WATER

www.EKULF.com