

ONE SHADING MEETS BULK FILLING



Venus Bulk Flow ONE.

The first flowable one-shade bulk-filling composite for everyday cases.

Always grab the right shade: With the incredible shade-matching properties of Venus Bulk Flow ONE the restoration becomes invisible by blending into the surrounding dentition.

Efficient handling: Cavities are filled quickly and easily with increments of up to 4mm without the need for an additional capping layer.

Long-lasting restorations: Low shrinkage stress, high flexural strength and wear resistance ensure durable restorations that withstand the mastication load in the posterior region.



Try it out for yourself: kulzer.com/try-venus-bulk-flow-one



KULZER
MITSUI CHEMICALS GROUP

4
3
2
1
0
mm

ONE LAYER FOR ALL SHADES




Venus Bulk Flow ONE simplifies and streamlines the restorative process.

Venus Bulk Flow ONE combines the self-levelling effect of a flowable material with a high curing depth and good stability after application. This allows to place increments of up to 4 mm inside a cavity.


With Venus Bulk Flow ONE dentists can reduce treatment time without the need for shade selection or the placement of a universal composite as an additional capping layer.

With its range of indications from bulk filling through cavity lining to fillings in primary teeth, the efficient restoration procedure is further enhanced by the equally reliable and versatile bonding solution iBOND Universal.

VERSATILITY IN APPLICATIONS

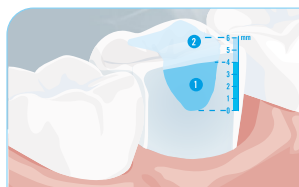

BULK & BODY TECHNIQUE:
 1 Venus Bulk Flow ONE
 2 + e.g. Venus Diamond/Pearl

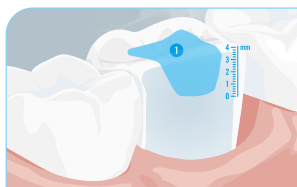

BULK FILLING TECHNIQUE:
 1 Venus Bulk Flow ONE


EFFICIENT HANDLING



Discoloured cavity



Venus Bulk Flow ONE
after polishing

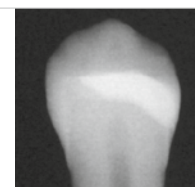

High radiopacity
ensuring a reliable
diagnosis

Art. code	Package	Content
Venus Bulk Flow ONE – Refills		
66094723	VENUS BULK FLOW SYR 1x2G ONE	1x 2g syringe 5x Cannulas
66094724	VENUS BULK FLOW PLT 2x 10x0.20G ONE	2x 10x0.20g PLT
Venus Bulk Flow ONE – Value Packs		
66095335	VENUS BULK FLOW ONE VALUE-KIT SYRINGE	10x2g syringes 50 Cannulas Pictorial card
66095388	VENUS BULK FLOW ONE VALUE-KIT PLT	80x0.20g PLT Pictorial card



Discover our new ONE: kulzer.com/try-venus-bulk-flow-one