

Pathfinding Files

-For Advanced Treatment-

MANI®

2 Hand Files For Pathfinding

D Finders



Best for calcified root canals


The D-shaped cross section of the D Finders increase blade strength, and prevent excessive filing and binding in calcified or curved root canals. Pre-curve the tip of the blade before use as a first step in the root canal negotiation.

Product name: MANI® D FINDERS

Glide Finders

For negotiating & cutting


The Glide Finders feature a varying cross section with multiple-taper blade, providing optimal tip resilience and shank flexibility near the handle.

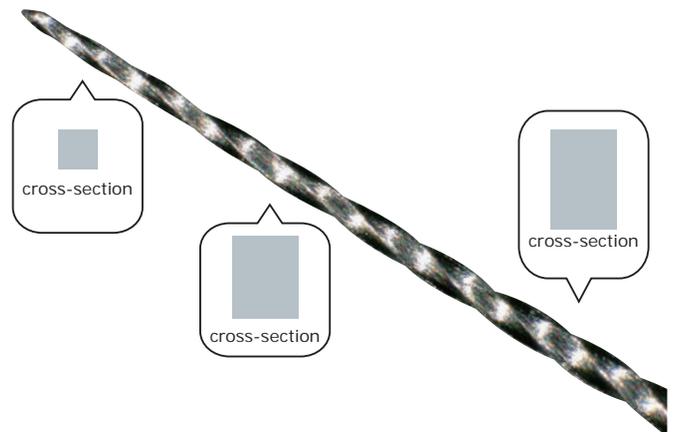

Product name: MANI® K-FILES

Both Files have high strength and are recommended when it is not possible to penetrate root canals by conventional K-FILES.

Size	Tip Diameter (mm)	Color code
08	0.08	
10	0.10	
12	0.12	
15	0.15	

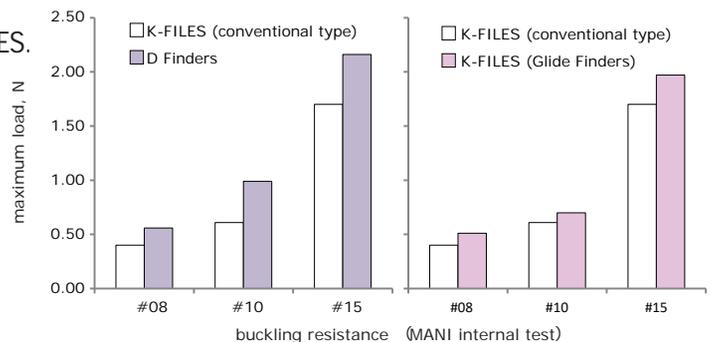

⚠ Please read IFU before use.

www.mani.co.jp

MANI, INC. 8-3 KIYOHARA INDUSTRIAL PARK, UTSUNOMIYA, TOCHIGI, 321-3231 JAPAN

Tel: +81 28 667 8592 E-mail: dental.exp@ms.mani.co.jp
URL: <http://www.mani.co.jp>



EC REP Obelis s.a.

Bd Général Wahis 53 1030 Brussels Belgium D10810